

Accurate Quantification with ICPL: The Quantification Study of the ABRF Proteomics Research Group 2006

The objective of the ABRF PRG study was to quantify 8 proteins in 2 mixtures. An ultraflex IITM was used to analyze ICPL-labeled protein digests using an LC-MALDI approach under control of the WARP-LCTM software. An average quantification error of 9.3 % was achieved using high performance instrumentation, Prespotted AnchorChips, ICPL technology and sophisticated bioinformatics tools - being the best result amongst all iTRAQTM and ICPLTM data returned to the study.

Sample

Eight proteins were supplied by ABRF PRG in 2 mixtures. BSA was specified as being present at 1:1 ratio. 3 unspecified proteins were also 1:1 “regulated”, for the other 4 proteins a maximum regulation of 1:100 was possible. Absolute protein amounts and ratios of the mixture were provided by the organizer after all experimental results had been submitted, these are shown in Figure 1.

Methods

ICPL-labeling of the 8 protein mixtures

The two protein samples were labeled with the ICPL (light or heavy, L or H) reagent after reduction and carbamidomethylation, as specified in the manual of the SERVA-ICPL-KitTM, using a protein concentration of 4 mg/ml. The concentrations of the individual proteins in each mixture varied from 3-600 pmol. After ICPL-

labeling, the two samples were combined and subjected to a proteolytic digest overnight using either trypsin (Promega) or endoproteinase Glu-C (Roche) following standard protocols. The digestion was stopped with 0.1% TFA.

LC and fractionation of samples

The Agilent CapLC 1100 was used for chromatographic separation controlled by the HyStar software. An aliquot of 8 μ l ICPL-labeled protein digest mixture was injected onto a C18-PepMap column (180 μ m id.x15cm; LC Packings) and separated with an 80 min gradient (5-35% acetonitrile in 0.1% TFA).

The eluting peptides were deposited directly onto a 384 Prespotted AnchorChipTM (PAC384) using the PROTEINEER fcTM fraction collector. The PAC targets come ready to use with prespotted matrix anchors, therefore, an additional matrix make-up flow was not required, and calibration standards are also pre-spotted, for easy external calibration. Prior to the MS-analysis the dried samples were washed for approx. 5 sec with 10mM NH₄H₂PO₄ (Sigma-Aldrich) in 0.1% TFA.

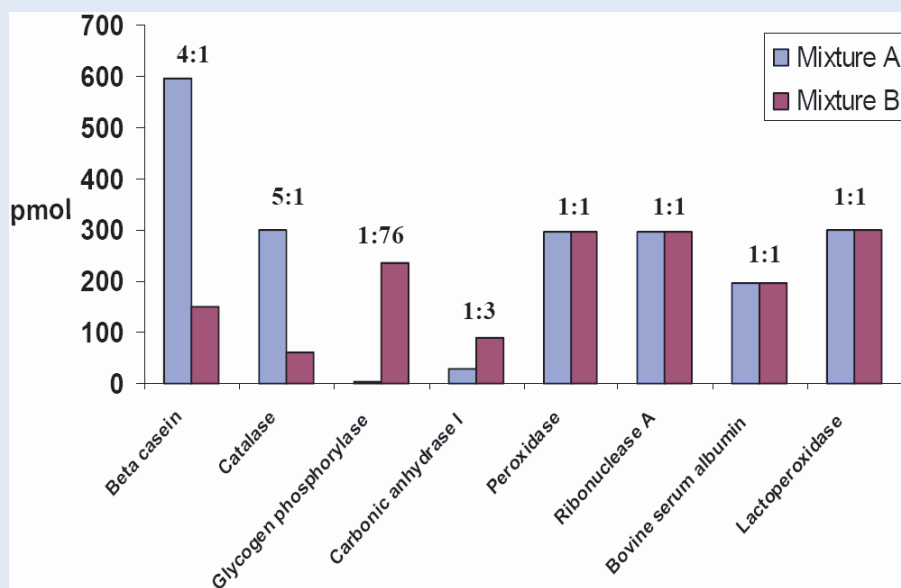


Fig. 1: Composition of the samples A and B analyzed in the PRG study.

Table 1: Result compilation from three LC runs as provided to the PRG study. Samples A and B were labeled with the ICPL forms L and H, respectively. CV: coefficient of variation between the individual protein quantifications in the three runs, which is typically less than 5 %; avg. L/H: experimental A/B ratios for identified proteins listed, Correct: theoretical L/H ratios; Error: relative error for the H/L value given for each protein.

| | | H/L | H/L | H/L | CV [%] | avg L/H | Correct | Error [%] |
|--|--|--------------|--------------|----------------|---------------------------|---------|---------|-------------|
| 1 | (P00432) Catalase | 0,17 | 0,18 | 0,17 | 2,7 | 5,67 | 5 | 13,4 |
| 2 | (P02663) Alpha-S2 casein precursor | 0,22 | 0,21 | | 1,6 | 4,62 | 4 | 15,5 |
| 3 | (P02666) Beta casein precursor | 0,22 | 0,23 | | 0,7 | 4,43 | 4 | 10,9 |
| 4 | Alpha-S1 casein precursor | 0,26 | 0,22 | | 8,1 | 4,16 | 4 | 4,0 |
| 5 | (P02668) Kappa casein precursor | | 0,23 | 0,26 | 5,7 | 4,05 | 4 | 1,2 |
| 6 | Ig heavy chain precursor (B/MT.4A.17.H5.A5) - bovine | 0,25 | 0,25 | | 0,8 | 4,03 | 4 | 0,8 |
| 7 | (P24627) Lactotransferrin precursor (Lactoferrin) | 0,98 | | 0,99 | 0,5 | 1,02 | 1 | 1,5 |
| 8 | (P80025) Lactoperoxidase precursor (LPO) | 1,01 | 0,99 | 1,03 | 1,6 | 0,99 | 1 | 1,0 |
| 9 | (P61824) Ribonuclease pancreatic (RNase A) | 1,08 | 1,01 | 0,98 | 4,1 | 0,98 | 1 | 2,3 |
| 10 | (P02769) Serum albumin precursor (BSA) | 1,05 | 1,06 | 1,00 | 2,5 | 0,96 | 1 | 3,5 |
| 11 | (P00915) Carbonic anhydrase I (CA-I) | 3,22 | 3,32 | 3,52 | 3,8 | 0,30 | 0,33 | 9,7 |
| 12 | (P00489) Glycogen phosphorylase, muscle form | >50 | >50 | 96,88 | | 0,010 | 0,013 | 21,6 |
| 3 LC-MALDI runs from 2 digests: | | Glu-C | Glu-C | Trypsin | Average Error [%]: | | | 7,11 |

MS and MS/MS data acquisition

All spectra were obtained fully automatically on a MALDI-TOF/TOF (ultraflex II TOF/TOF™). MS spectra were acquired at constant laser power and externally calibrated with the peptide calibrant mixture provided on the PACs. After the MS-analysis of the complete LC-MALDI run, the software WARP-LC defined a list of unique compounds and selected best suited precursor ions for MS/MS data acquisition. All MS/MS spectra were acquired for each chromatographic compound satisfying quality criteria (S/N etc.). H/L ratio thresholds were not applied as criteria for selecting MS/MS precursors in order to quantify 1:1 regulated proteins. In real experiments such regulation thresh-

holds would be used to significantly reduce analysis time for LC-MALDI-MS/MS. The Glu-C digest was analyzed twice by LC-MALDI and the trypsin digest once. After data acquisition of each LC-run a MASCOT database search was triggered.

Quantification and data validation

All peaks were picked using the proprietary SNAP algorithm that fits isotopic patterns to the matching experimental data even under conditions of significant overlap of isotopic peak clusters. The algorithm provides the monoisotopic mass, the intensity and area under the envelope of the isotopic cluster and the resolution of the peaks in the cluster. In particular, peaks of lower intensity in general suffer noise interferences and, there-

fore, reduced reliability of the determination of peak resolution. For best quantification results under conditions of high dynamic range (1:100 was the given dynamic range in this study), peak intensities and not areas were used for quantification. After acquisition of MS and MS/MS spectra and database searching, fully automatic quantitative analyses were calculated during the Mascot result compilation in the WARP-LC ProteinBrowser™ (see Fig. 4). The average and median H/L values plus the CV values were the major parameters used for the validation of the data. Additional parameters that were included in the interactive result acceptance process were the relative intensity and the Mascot scores of the individual peptides.

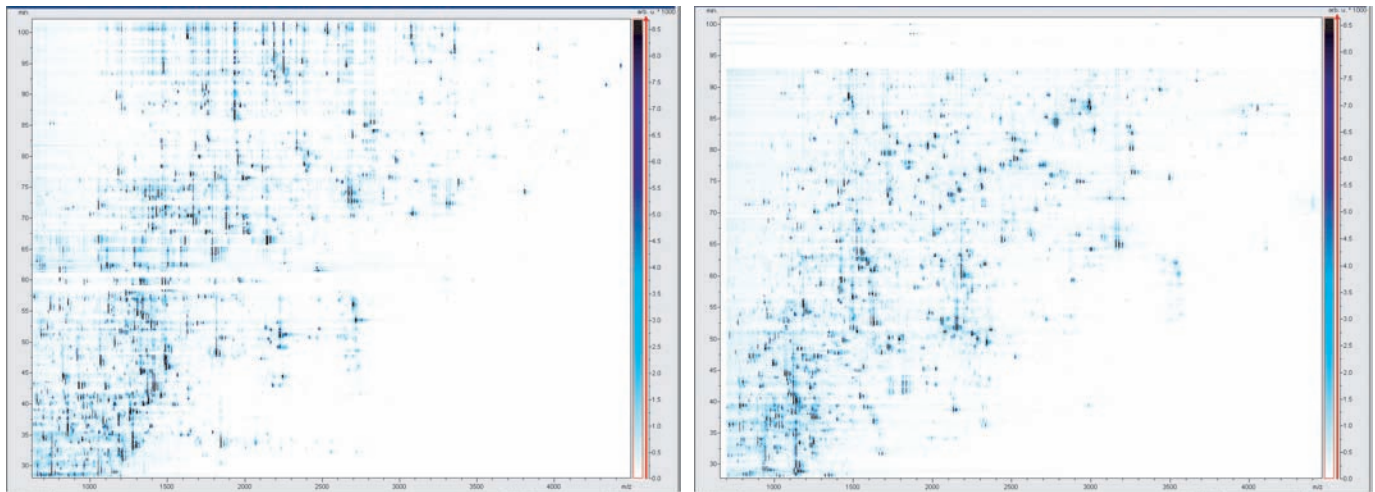


Fig. 2: Survey Views of the entire LC-MS runs of the Trypsin (left) and the Glu-C (right) digests of the ICPL labeled samples showing similar numbers of peptides and similar molecular weight distribution.

Table 2: Detailed compilation from the ABRF PRG study for participants using commercially available quantitation chemistries (iTRAQ, ICPL), sorted by increased average % error across all 8 proteins. Code: anonymous identifier of participants – v indicates participating vendors, 270960v indicates the results of Bruker Daltonics (BDAL); Label: used chemistry; Instrument: used MS technology; Mode: quantification based on area or peak height. Colors are used to guide the eye to similar groups in terms of performance or methods.

| Code | Label | Instrument | ID Software | Quant Software | Mode | Average Error [%] | |
|--------|---------------|------------|----------------|----------------------|--------------------------|-------------------|------------|
| 27960v | BDAL | ICPL | TOF/TOF | Mascot | WARP-LC 1.1 | Peak Height | 9,3 |
| 98765 | | iTRAQ | q-TOF; TOF/TOF | Mascot | Excel | Peak Height | 19 |
| 11576 | | ICPL | TOF/TOF | Mascot | GPS Explorer 3.5 | Peak Area | 20 |
| 12125v | Vendor | iTRAQ | LTO | Sequest | developmental software | Peak Height | 25 |
| 21806 | | ICPL | TOF/TOF | Mascot | Excel; GPS Explorer | Peak Area | 54 |
| 55000 | | iTRAQ | QTRAP | ProID | ProQuant | | 79 |
| 51892 | | iTRAQ | q-TOF | Mascot | ProQuant; Excel | Peak Area | 85 |
| 15024 | | iTRAQ | TOF/TOF | Mascot | GPS Explorer 3.5 | Cluster area | 98 |
| 55001 | | iTRAQ | q-TOF | Mascot | Excel | | 105 |
| 12569 | | iTRAQ | q-TOF | Mascot | ProQuant | | 119 |
| 72911 | | iTRAQ | q-TOF; TOF/TOF | Mascot; Water's PLGS | Waters EAP; GPS Explorer | Peak Area | 142 |
| 65101v | Vendor | iTRAQ | TOF/TOF | Mascot | GPS Explorer | Peak Area | 164 |
| 29471 | | iTRAQ | TOF/TOF | Mascot; Phenyx | Excel | Peak Height | 195 |

Results

More proteins were obtained from the Glu-C digest as compared to the trypsin digest (see Tab. 1) although differences in terms of numbers and molecular weight distributions of the peptides were not obvious (Fig. 2). The sequence coverage of identified proteins was found to be dependent upon the individual protein rather than the protease with which it was digested. However, this result suggests that in any approach to quantification it is important to use more than one protease to ensure sufficient peptide pair numbers for reliable quantification

covering a greater number of proteins. Glu-C certainly is one attractive selection for that purpose. Across all identified proteins and including data from proteins which were contaminants of the expected proteins, the absolute average error of quantification was found to be 7.1 % (Tab. 1).

Several proteins were identified in addition to the 8 proteins included in this study, which were obviously contaminants of the proteins from the commercial source and which provided identical quantitative readings as the intended specified proteins: e.g. β -casein (L/H = 4) is accompanied by

additional α -S1 and α -S2 and β -casein found with identical L/H value.

A more detailed comparison of the quantification results was obtained using data compiled by the organizers of the PRG study. This comparison examined the data quality obtained from different labeling strategies using standard commercial kits, i.e., ICPL-Kit (Bruker Daltonics/SERVA) and iTRAQ (ABI/Invitrogen) (Tab. 2). Some participants also used C-terminal $^{18}\text{O}/^{16}\text{O}$ labeling. Interestingly, ICATTM or iPROTTM were not used at all.

Table 3: Detailed result compilation from Table 2, sorted by increased error level across all 8 proteins. All tabulated values are absolute % errors of quantitation of individual proteins as well as average values across all proteins (Dynamic range <100) or excluding Glycogen phosphorylase (Dynamic range < 1:10). Error "0" typically means that participants identified a 1:1 regulation and then set the reported ratio deliberately to 1. Colors are used to guide the eye to similar groups in terms of performance.

| Code | BSA | Catalase | GlyPhosb | CA | Perox | Rnase A | LacPerox | β Cas | # of Proteins | Average (Dynamic range) | |
|--------|-----|----------|----------|-----|-------|---------|----------|-------------|---------------|-------------------------|--------|
| | | | | | | | | | | < 1:100 | < 1:10 |
| 27960v | 4 | 13 | 24 | 10 | | 2 | 1 | 11 | 7 | 9,3 | 5,9 |
| 98765 | 0 | 20 | | | 0 | 0 | 0 | 97 | 6 | 19 | 17 |
| 11576 | 3 | 1 | 128 | 1 | | 0 | 6 | 3 | 7 | 20 | 2,0 |
| 12125v | 26 | 33 | 52 | 10 | 23 | 17 | 23 | 17 | 8 | 25 | 19 |
| 21806 | 0 | 3 | 280 | 61 | | 6 | 7 | 22 | 7 | 54 | 14 |
| 55000 | 3 | 30 | 432 | 68 | | 2 | 1 | 19 | 7 | 79 | 18 |
| 51892 | 0 | 85 | 24 | 59 | 23 | 400 | 20 | 73 | 8 | 85 | 82 |
| 15024 | 5 | 23 | 584 | | 15 | 6 | 2 | 54 | 7 | 98 | 15 |
| 55001 | 13 | 21 | 584 | 38 | | 5 | 3 | 73 | 7 | 105 | 22 |
| 12569 | 7 | 17 | 660 | 1 | 56 | | 5 | 90 | 7 | 119 | 25 |
| 72911 | 0 | 60 | 660 | 200 | 0 | 70 | 96 | 50 | 8 | 142 | 60 |
| 65101v | 0 | 40 | 812 | 83 | | | 11 | 40 | 6 | 164 | 29 |
| 29471 | 0 | 37 | 1040 | | 8 | | 9 | 73 | 6 | 195 | 21 |

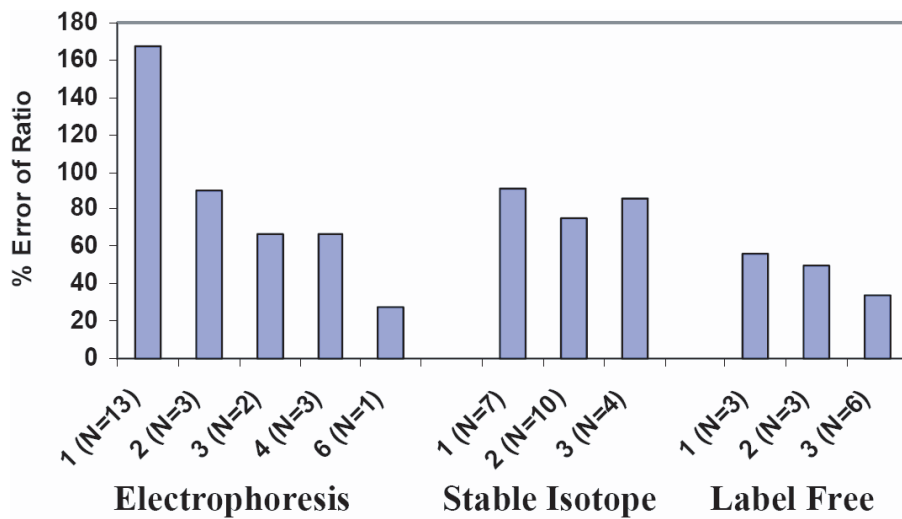


Fig. 3: Average % errors for the different analysis methods as a function of experiment repetition and quantification technique.

Typical average errors of the protein quantifications across all participants are listed in Table 2 in the red box. In Table 3 the average errors obtained for the individual proteins are listed. The errors range from 5 % for BSA (this was the known 1:1 standard) through 20-50 for most proteins to > 400 % for the highly regulated Glycogen Phosphorylase, H/L = 1:76. Across all proteins the average quantification results for the individual participants ranged from 9.3 to 195 % error margins and from 8 to 6 identified proteins (only results > 5 IDed proteins are

listed here).

The Bruker Daltonics result using ICPL and dedicated quantification software (WARP-LC 1.1) is coded 270960v and ranked at the top of the quality sorted list. In addition, the ICPL chemistry was used by 3 participants within the top 5, indicating its usefulness for quantitative proteomics. Surprisingly, the iTRAQ-labeling approach using the standard quantification software that was written for it, provided dramatically high average quantitation errors of 80-200%, even

Microsoft Excel based interpretation or developmental software provided superior results. Another interesting aspect was that in the top performing group of participants, the quantification based on the peak height rather than peak area was preferred.

Compared to the average results from all other quantification strategies, our results of 9.3 % average error was very competitive (see Fig. 3). We identified proteins with H/L = 1 with errors < 5 % and proteins with 2 < H/L < 10 with approx. 10 % errors. The quality of the quantification results was not significantly increased by repeated LC-runs (see Tab. 1) in agreement with Figure 3. This indicates a high level of confidence of the individual quantification results (CV values of quantification typically less than 5 % between LC runs, Tab. 1), provided sufficient numbers of peptides for calculation. This robustness is required if more complex samples on true proteomics scale are to be analyzed. On this scale the suitability of the ICPL chemistry for protein pre-fractionation in case of high proteome complexity will become important as high repetition numbers (>>2) need to be avoided due to outrageous workload problems.

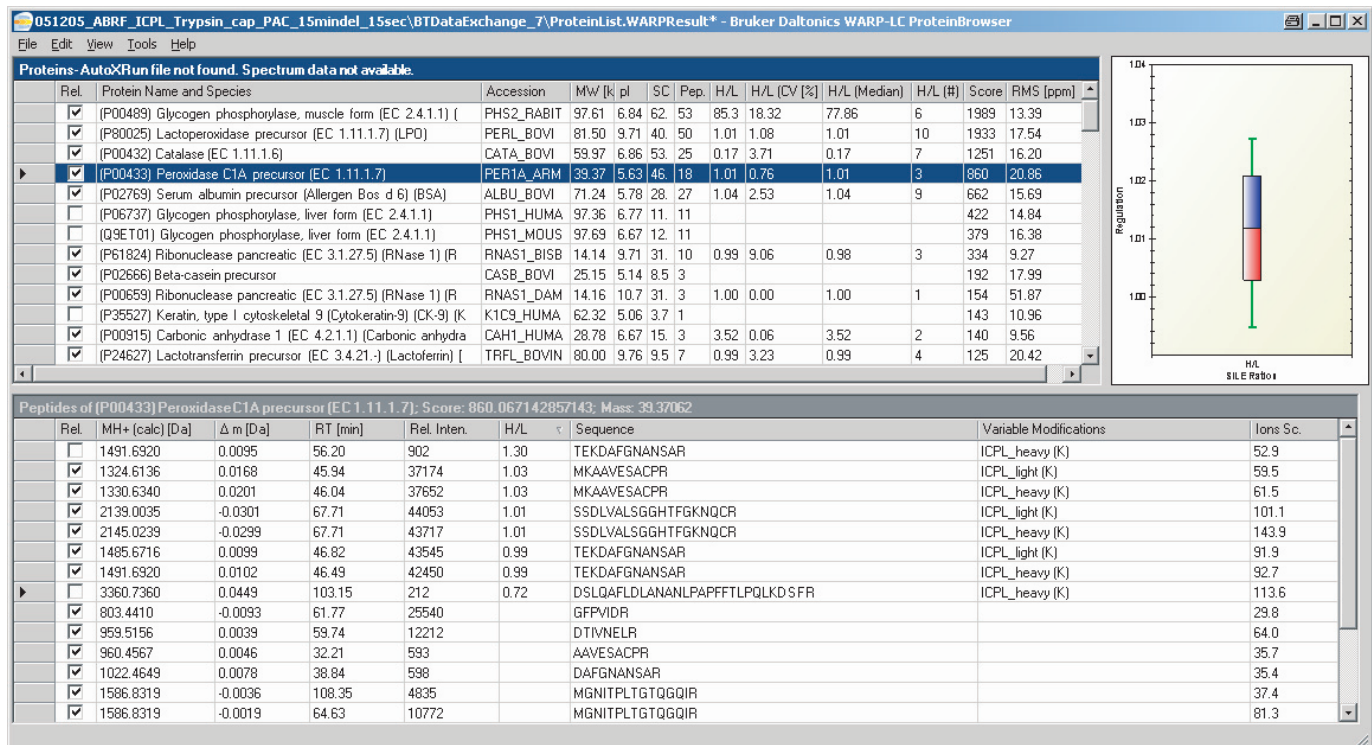


Fig. 4: WARP-LC ProteinBrowser displays the identification and quantification results and allows the manual acceptance of the results or validation on the level of peptides and proteins. Peroxidase has been selected from the protein table, the corresponding peptides and regulations are shown in the peptide list. Horseradish peroxidase was IDed and quantified on the 1 % CV level.

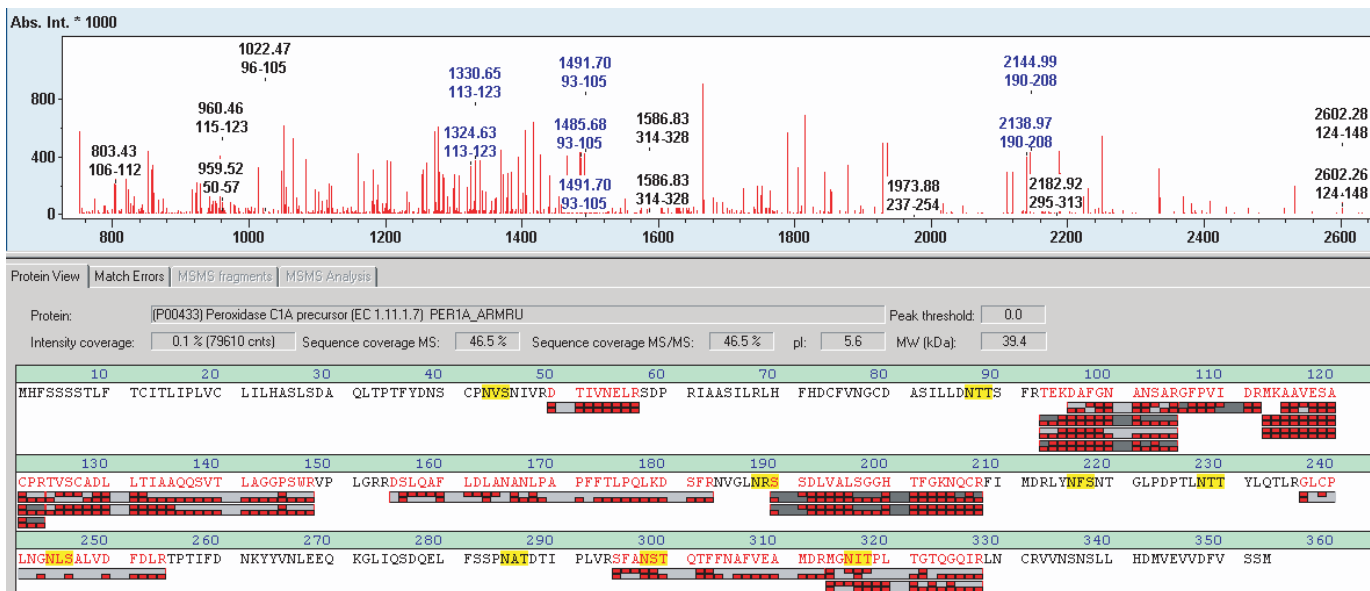


Fig. 5: BioTools Sequence coverage map of horseradish peroxidase. Matching peptides are shown as grey bars, darkness indicating peak intensity. The red bricks indicate the matching N-terminal (upper row) and C-terminal (lower row) fragment ions. Putative glycosylation sites are coded in yellow.

Data Analysis

We used the advanced data validation tools of WARP-LC 1.1 to determine reliable quantitative data (Fig. 4). For each identified protein a number of protein chemical parameters are

listed in the visualization tool (WARP-LC ProteinBrowser) of the software including the total peptide score, protein MW and pI, sequence coverage (%) and the number of IDED peptides. Information relevant to quantification includes the average H/L ratio, the number of peptide pairs utilized

for quantification, the coefficient of variation (CV: the ratio of the standard deviation to the mean) and the median of all values for a particular protein. Validation is typically based on the CV plus the difference between the average H/L value and the median. The box-and-whiskers plot is

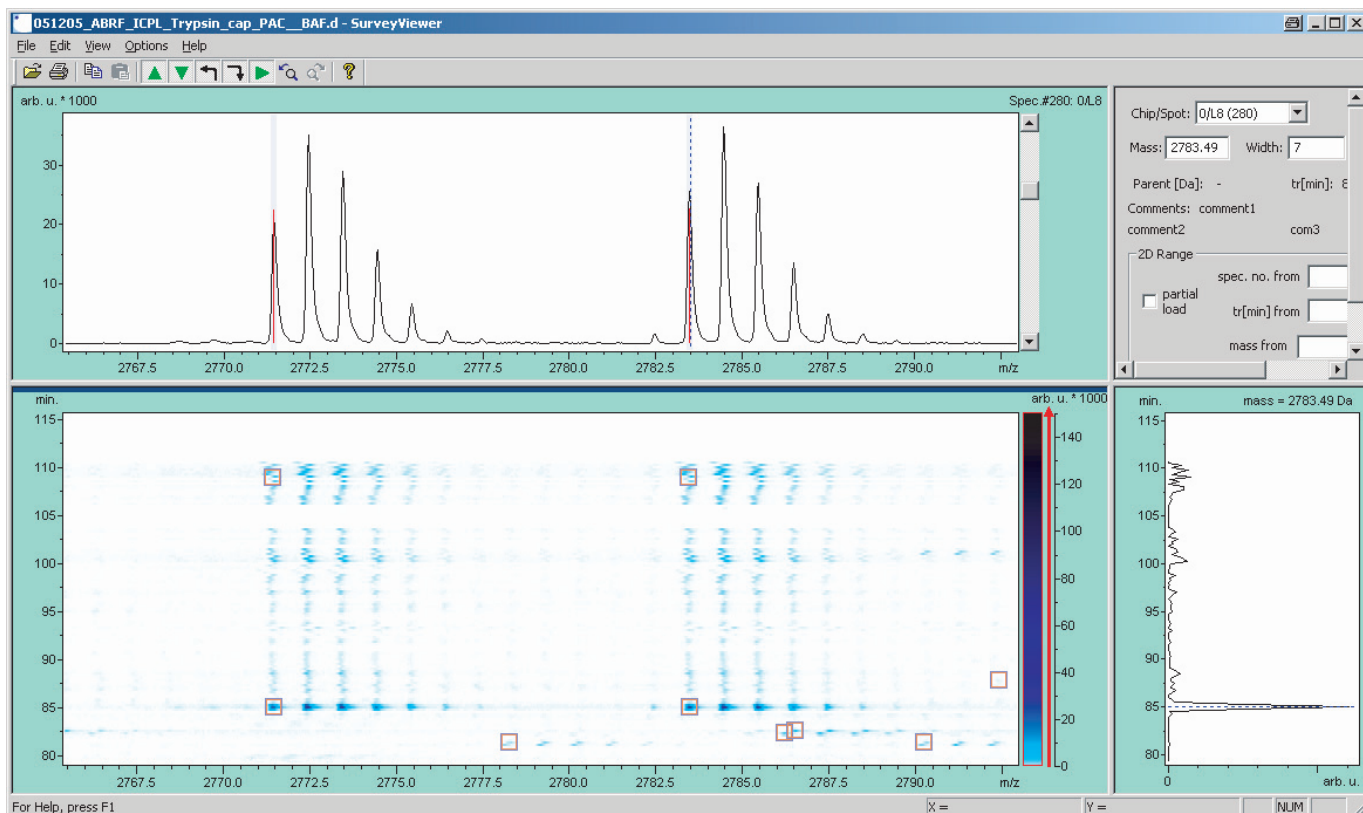


Fig. 6: Zoom View into the LC-run of the trypsin digest (see Fig. 2) that displays the 1:1 regulated pair of ions from lactoperoxidase. All selected precursors for MS/MS are labeled with squares by the WARP-LC software; the MS spectrum and extracted ion chromatogram for the 85 min fraction with m/z 2783.49 are shown on top and on the right side of the 2D display, respectively. MS/MS spectra of m/z 2783.49 shown in Fig. 7 were selected at 85 and 108 min.

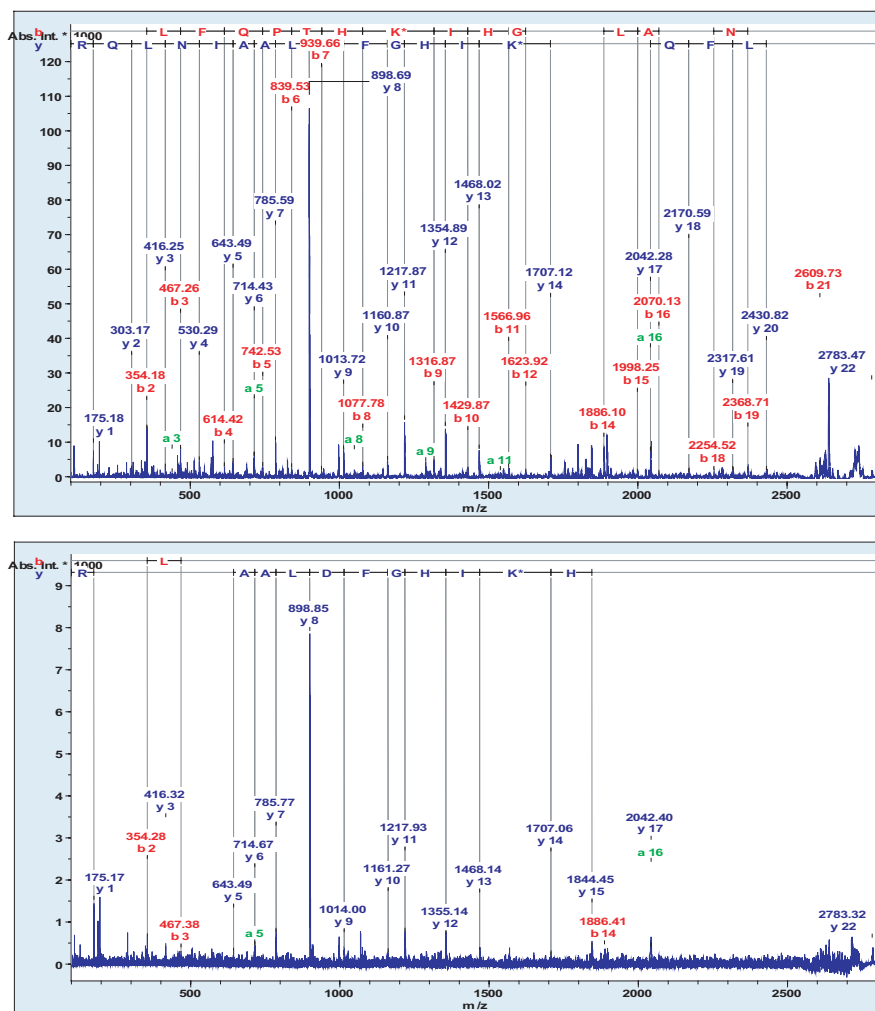


Fig. 7: BioTools 3.0 display of the m/z 2783.49 peptide from lactoperoxidase at 85 min (top) and 108 min (bottom) that were both identified automatically in the course of the analysis.

also displayed, which visually describes the maximal and minimal data point that is not detected as an outlier, the median and the body of the central 50 % of the data, i.e., the box.

In Figure 4, the WARP-LC Protein-Browser view highlights the peroxidase data, where median as well as the average H/L ratio were calculated as $H/L = 1.01$ (1.0 expected) with a CV of 0.8 %. Typical precursor mass accuracies achieved across all detected peptides of the identified proteins were in the 5-10 ppm range. Although horseradish peroxidase was found (P00433) (Figs. 4 and 5) it was not provided in the result list due to a transmission error to the ABRF website during the results upload via the web.

The ability to manually judge the quantification results and to select or deselect individual peptides to include or exclude them for the calculation of average/median values of protein

regulation (see Fig. 4), was an incredibly valuable tool for the high quality data provided in this study.

The high analytical sensitivity of LC-MALDI MS/MS was shown by the peptide m/z 2783.49 (matching lactoperoxidase) detected at 85 min then later at 108 min in a small maximum of a peak smear (Fig. 6). In the WARP-LC software, the automatic picking of precursor ions in such smears can be controlled. In the example, reacquisition of the same peptides was disabled for a long time, however after 25 min and after interruption of the ion trace along the chromatogram, the same parent mass was readily accessible for further MS/MS analysis even though the latter peak was represented by a weak ion. At both times, the same peptide was identified and the matching fragments are displayed (Fig. 7).

Conclusions

The analytical platform chosen for this study consisted of a number of unique components that laid the foundation for a great performance in the PRG-ABRF '06 quantification study (<http://www.abrf.org/index.cfm/group.show/Proteomics.34.htm>). An average error level $< 10\%$ at a dynamic range of up to 1:100 was the outstanding result of all stable isotope labeling results reported in the study and far beyond the average performance of gel, labeling or label free approaches reported in this study. The components of this analytical platform were:

1. **A high performance MALDI-TOF** with sub fmol MS/MS sensitivity and a high performance precursor ion selector (> 750 resolving power). The spectra acquisition under fuzzy logic control prevented peak saturation and large quantification errors for strongly regulated peptides and proteins.
2. **Prespotted AnchorChips** as single-use MALDI targets with zero memory effect (peptide background from previous LC-runs) that could interfere with high quality quantification and a high reproducibility of the ion current from different spots within one chromatographic fraction.
3. The **ICPL label technology** providing isotopically pure labels that work well even for proteins that are greatly regulated. In contrast to several other methods used in the study (iTRAQ, ion current), ICPL is fully compatible with quantification on the proteomics level in conjunction with protein pre-fractionation.
4. **Bioinformatics** with numeric as well as visual and interactive validation tools that are unique in the market and that ensure high quality results as the scientist has full control over the analytical process without the need to work with external tools such as spreadsheet calculation.

References

[1] PRG-ABRF '06 quantification study (<http://www.abrf.org/index.cfm/group.show/Proteomics.34.htm>)

Authors:

*Detlev Suckau, Stephanie Hahner,
Laura C. Main, Wolfgang Jabs*

*1: Bruker Daltonik GmbH, Bremen,
Germany*

*2: Bruker Daltonics Ltd., Coventry,
United Kingdom*

Keywords:

Quantitation of proteomes, ICPL,
Stable Isotope Labelling (SILE),
MALDI-TOF/TOF

Instrumentation & Software:

Agilent CapLC 1100 (Agilent),
C18-PepMap column (LC Packings),
HyStar, SERVA ICPL-Kit,
Prespotted AnchorChip (PAC384),
PROTEINEER fc,
ultraflex II TOF/TOF, WARP-LC,
BioTools (all Bruker Daltonics)

The SERVA ICPL-Kit is available from BDAL's CARE online portal (www.bdal.de/care).

| Product | Cat. No. | Size |
|-----------------|----------|-----------------------|
| SERVA ICPL™-Kit | #234017 | 1 Kit (2 x 6 samples) |

WARP-LC, HyStar, Prespotted AnchorChip, PROTEINEER fc, ultraflex II TOF/TOF and BioTools are trademarks of Bruker Daltonics corporation. ICPL is a trademark of TopLab GmbH, Germany. iTRAQ is a trademark of Applied Biosystems. ICAT is a trademark of the University of Washington. iPROT is a trademark of Agilix Corporation Inc.

www.bdal.de
www.bdal.com

BRUKER
DALTONICS®

Bruker Daltonik GmbH

Bremen · Germany

Phone +49 (421) 22 05-0
Fax +49 (421) 22 05-104
E-Mail sales@bdal.de

Bruker Daltonics Inc.

Billerica, MA · USA

Phone +1 (978) 663-36 60
Fax +1 (978) 667-59 93
E-Mail ms-sales@bdal.com

See our website for worldwide sales and service contacts.